Littmann® CORE Stethoscope

Connects to Eko Software

(en)	English	pp. 1	(sl)	Slovenski	pp. 33
(fr)	Français	pp. 3	(et)	Eesti	pp. 35
(de)	Deutsch	pp. 5	(lv)	Latviesu	pp. 37
(it)	Italiano	pp. 7	(lt)	Lietuviskai	pp. 39
(es)	Español	pp. 9	(ro)	Romana	pp. 41
(nl)	Nederlands	pp. 11	(hr)	Hrvatski	pp. 43
(sv)	Svenska	pp. 13	(bg)	Български	pp. 45
(da)	Dansk	pp. 15			
(no)	Norsk	pp. 17			
(fi)	Suomi	pp. 19			
(pt)	Português	pp. 21			
(el)	Ελληνικά	pp. 23			
(pl)	Polski	pp. 25			
(hu)	Magyar	pp. 27			
(cs)	Čeština	pp. 29			
(sk)	Slovensky	pp. 31			

3M™ Littmann® CORE Digital Stethoscope

Quick Start Guide

Charge the Littmann CORE Digital stethoscope before you first use it. After it's charged, it's ready to use.

eko

Download the Eko App (where available) from the Apple App Store® or the Google Play Store® to access additional features* of the Littmann CORE stethoscope including wireless listening, advanced sound and filtering control, murmur detection, saving recordings, and sharing PDF reports.

*Select features may require a paid subscription and may not be available in all regions or all languages.


Scan with your phone's camera to download the Eko App.

In the Box

- Eko CORE Digital Attachment and 3M™ Littmann® Cardiology IV™ Stethoscope
- Additional set of 3M Littmann Snap Tight Soft-Sealing Eartips (small)
- Nonchill Bell Sleeve
- USB Charging Cable
- Quick Start Guide

Eartips

Insert eartips in a forward direction.

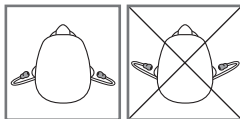Volume

Use the (+) and (-) buttons to adjust volume. There are 7 volume steps.

Turn On and Off

When the Littmann CORE Digital stethoscope is off, it functions as an unamplified stethoscope.

To turn the Littmann CORE Digital stethoscope on, press the power toggle in. The light will begin flashing. You will now hear amplified and filtered sound in the eartips. The Littmann CORE Digital stethoscope must be on to use the digital features.

Active Noise Cancelling

The Littmann CORE Digital stethoscope includes active noise cancelling to remove ambient noise from the surrounding environment. Active noise cancelling can be turned off or on through the Eko App.

Updates and New Features

Eko releases new features to improve your stethoscope. Connect it to the Eko App to get the latest device updates.

Telemedicine and Analysis Features

The Littmann CORE Digital stethoscope is compatible with premium telemedicine and analysis features (where available, sold separately) offered by Eko. Learn more at ekohealth.com.

Notifications

(blinking white)

CORE is on & seeking device

(solid white)

CORE is on & connected

(solid blue)

CORE is playing back a recording

(blinking orange)

CORE is low on battery

(solid orange)

CORE is off & charging

(solid green)

CORE is fully charged

Recording

You can use the Eko App to take recordings when connected to the Littmann CORE Digital stethoscope. To learn more, go to support.ekohealth.com.

Charging

Connect the micro USB cable to the Littmann CORE Digital stethoscope and use a UL-certified or CE-marked 5V USB wall charger. The Littmann CORE Digital stethoscope will take about 2 hours to fully charge. The light will turn green when it is fully charged.

User Manual

Review the user manual before using your Littmann CORE Digital Stethoscope. To view the manual visit: ekohealth.com/ifu

For more information on 3M devices, see HCBGRegulatory.3M.com

Service and Support

In the USA, please visit littmann.com/service or call 1-800-228-3957. If you are outside of the USA, please visit littmann.com for the contact information of your local 3M office.

Apple and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc., registered in the U.S. and other countries.

Google Play and the Google Play logo are trademarks of Google LLC.

The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by Eko is under license.

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System assembled by 3M Company

3M Deutschland GmbH
Health Care Business

3M (Schweiz) GmbH
Eggstrasse 91
8803 Rueschlikon
Switzerland

3M United Kingdom PLC
3M Centre, Cain Road
Bracknell, RG12 8HT,
United Kingdom

www.littmann.com

Stethoscope Made in USA with
Globally Sourced Materials

3M Company
2510 Conway Ave.
St. Paul, MN 55144 USA

3M Deutschland GmbH
Health Care Business
Carl-Schurz-Str. 1
41453 Neuss, Germany

3M (Schweiz) GmbH
Eggstrasse 91
8803 Rueschlikon
Switzerland

UK Responsible Person:
3M United Kingdom PLC
3M Centre, Cain Road
Bracknell, RG12 8HT,
United Kingdom

Eko CORE Digital Attachment
Made in Malaysia

Eko Health, Inc.
2100 Powell Street Suite 300
Emeryville, CA 94608 USA

Emergo Europe
Westervoortsedijk 60
6827 AT Arnhem
The Netherlands

MedEnvoy
Switzerland
Gotthardstrasse 28
6302 Zug
Switzerland

UK Responsible Person:
Emergo Consulting (UK) Limited
c/o Cr360 - UL International
Compass House,
Vision Park Histon
Cambridge CB24 9BZ
United Kingdom

©2023 3M; Eko Health, Inc. All rights reserved. 3M, Littmann, the L Littmann logo, and the L are marks and/or registered marks of 3M. Unauthorized use prohibited. Used under license in Canada. Eko trademarks, including CORE, are marks and/or registered marks of Eko Health Inc. Unauthorized use prohibited.

Issue Date: 2024-12

34-2001-0562-4 0191447.1

 Requester: Susan Barker Division Medical Solutions Division Agency project number: deZinnia_34844 PID: 16853 File Name: Insert_34-2001-0562-4_0191447-1.indd Material Spec ID: 34-2001-0562-4 Structure #: SS-85147 Date: 12/10/24	Printed Colors – Front: 
	Printed Colors – Back: 
	Match Colors: 
Scale:  1 Inch	This artwork has been created as requested by Solventum. Solventum is responsible for the artwork AS APPROVED and assumes full responsibility for its correctness.

PROOF REQUIRED

Please refer to Solventum Corporate General Spec. 205 for details